

merilas


shortpulse **810**

Laser therapy

for advanced and
end-stage glaucoma

your laser specialist

meridian 
medical

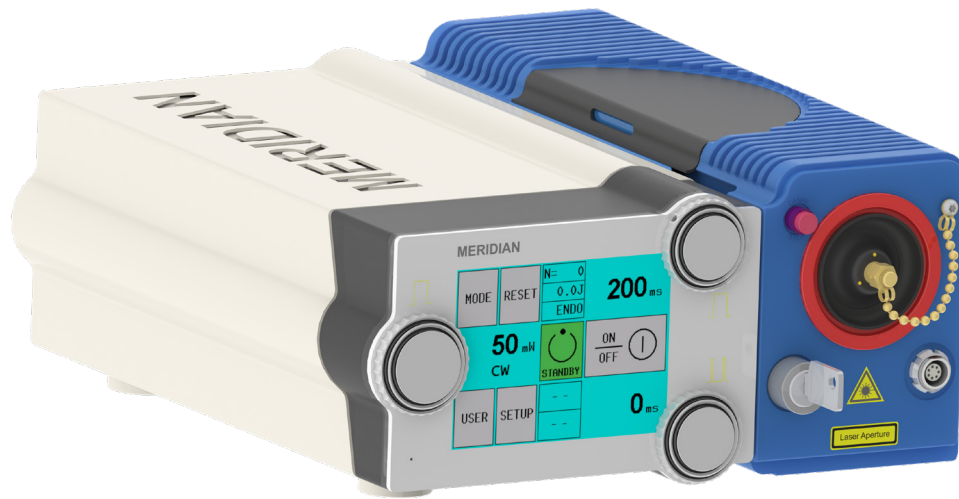shortpulse **810**shortpulse **810**

Diode Laser

Superior Excellence Quality Reliable Flexibility Comfortable Usability Material Medical Integrated Security Tested Longevity


your laser specialist

merilas
shortpulse **810**

SUPERIOR QUALITY & LONGEVITY

The Merilas 810 shortpulse is cooled by a unique thermal electrical cooling system. No ventilation slots are required. This ensures that no dust can penetrate into the interior of the laser.

The housing of the laser is made of a high grade aluminium giving the unique solid feeling of Meridian products while moving the laser and protects the sensitive parts of the laser.

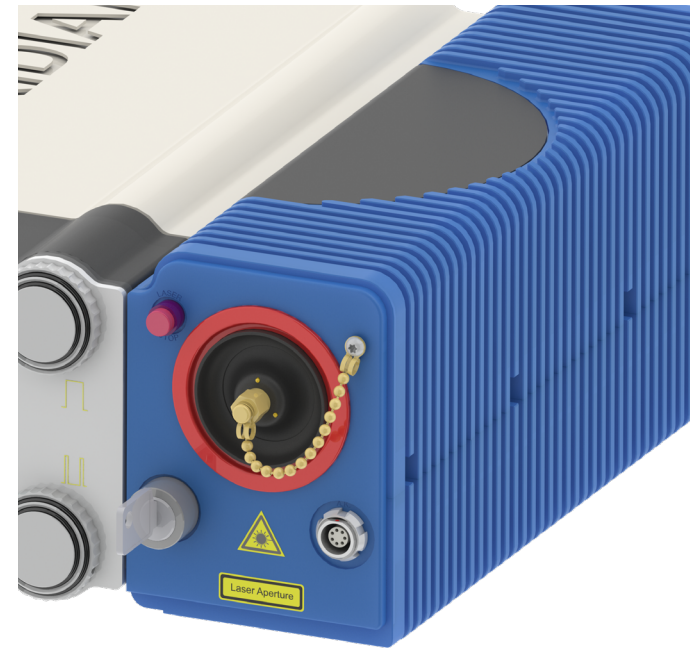

merilas


shortpulse **810**

USABILITY

The Merilas lasers convince by their fast start-up.

The detachable touch display with glass technology ensures flexibility and cleanliness at the workplace. The user interface is easy to use thanks to its intuitive design. Due to our special and unique thermal-electrical cooling system there are no disturbing noises or air turbulences.

All Merilas lasers impress with their small and compact size and are easy to transport. Each laser comes with a robust and practical transport case.

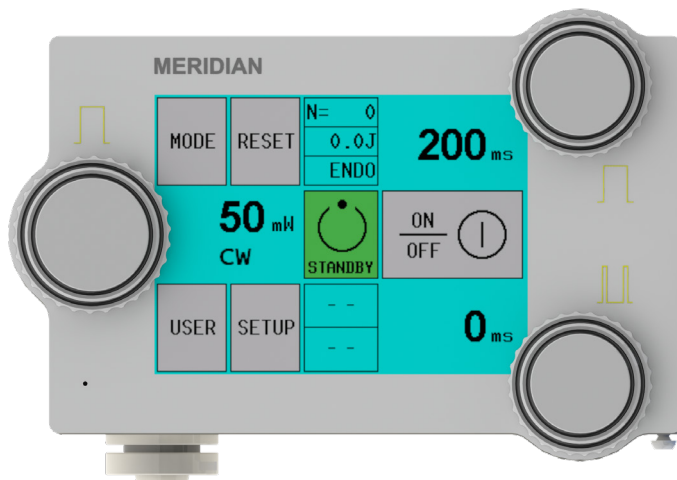

merilas

shortpulse **810**

SAFETY

The separate menu modes for easy identification of continuous wave and shortpulse mode ensure safe handling.

Our color-coded laser coupling to match the wavelength output makes it easy to assign the laser delivery systems. Auto Key connector: The automatic recognition of the laser delivery device must be confirmed by the user.

Our technicians can support you via remote service in case you need assistance. This function allows fast and safe troubleshooting.



merilas

shortpulse **810**



FLEXIBILITY & COMFORT

One laser: Two clinical applications (Glaucoma and Retina treatments).

Due to the small and compact design of the laser and its handy carrying case, the Merilas 810 shortpulse can easily be transported to other practices or clinics.

Our technicians can support you via remote service in case you need assistance. This function allows fast and professional troubleshooting.

The history of Meridian AG, now showing up as Meridian Medical Group, and the history of the medical Nd:YAG laser are closely connected. Meridian AG was already significantly involved in the development of the “Merilite” and 2006 the first Merilas laser of this laser family.

For the shortpulse lasers, new technology was developed and patented by our developmentengineers.

We select and integrate the best Swiss and European laser components to ensure the highest level of quality and long term reliability.

We use tested and reliable best practices in engineering and integration, ensuring highest performance in each of our systems.

Our highly skilled and experienced staff works to deliver the service and results our customers deserve and have come to expect.


- Not using your laser over a longer period of time will shorten its lifetime
- Store the laser in the supplied case if the laser is not used for a longer time
- Regular cleaning ensures stable operation of the laser



merilas

shortpulse **810**

CLINICAL INDICATION

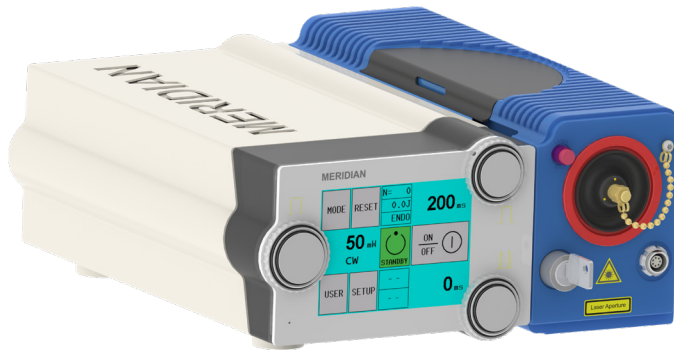

Retinal Photocoagulation:


- Retinal Photocoagulation with endo-probe
- Diabetic Retinopathy
- Central Serous Chorioretinopathy
- Diabetic Macular Edema
- Central/Branch Retinal Vein Occlusion
- Retinal Tears, Holes and Detachments

Glaucoma:

- Cyclophotocoagulation


WAVELENGTH BENEFITS – WHY 810 NM?

- The 810 nm wavelength is less absorbed by melanin and haemoglobin
- This wavelength has the advantage of penetrating the sclera and excessive haemorrhage

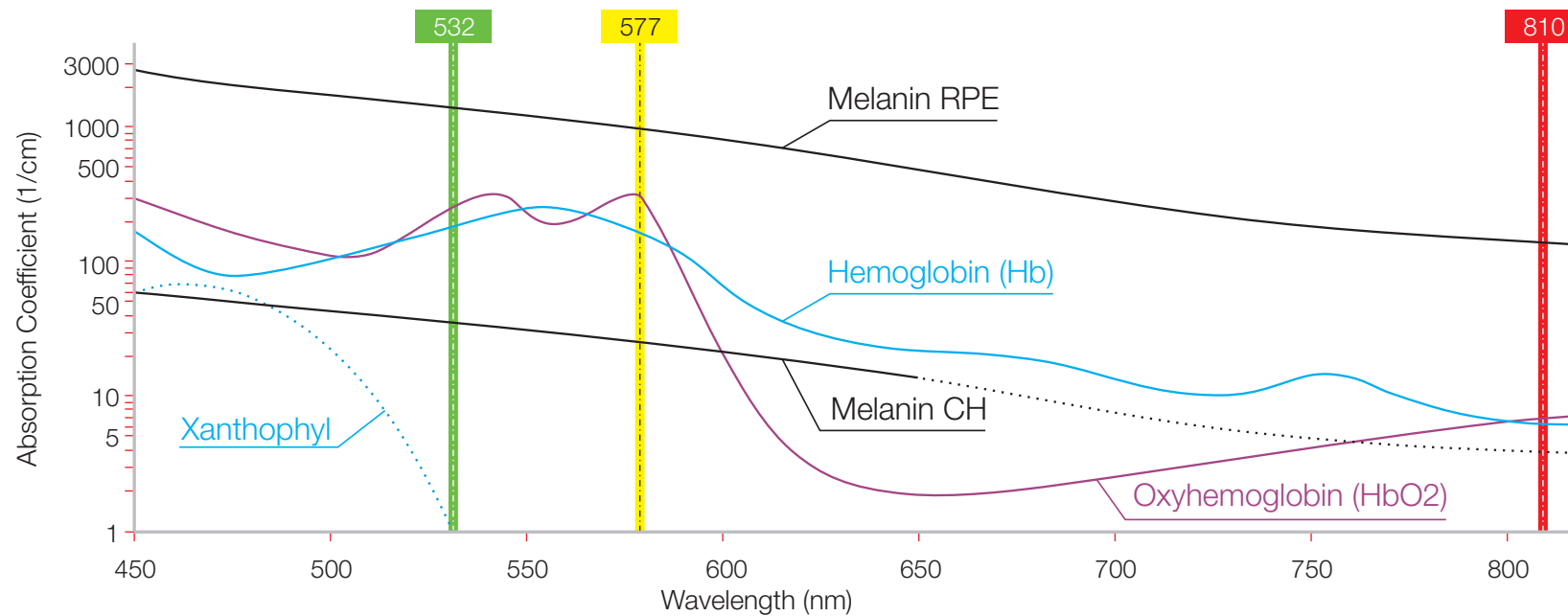


merilas

shortpulse **810**

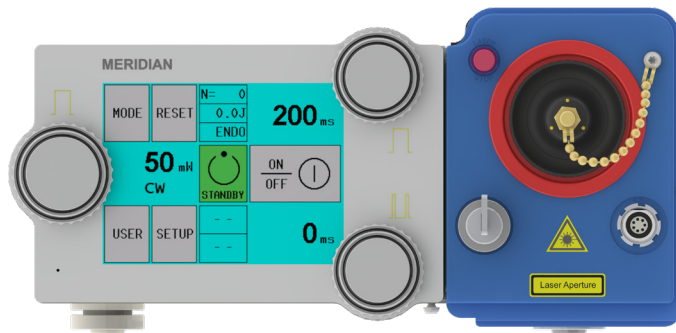

THE PRINCIPLE OF SHORTPULSE

- In shortpulse mode a pulse duration consists of many alternating short bursts and intervals.
- In contrast to the continuous wave (CW) mode, the tissue is not heated very much in the shortpulse mode – tissue is treated more gently.

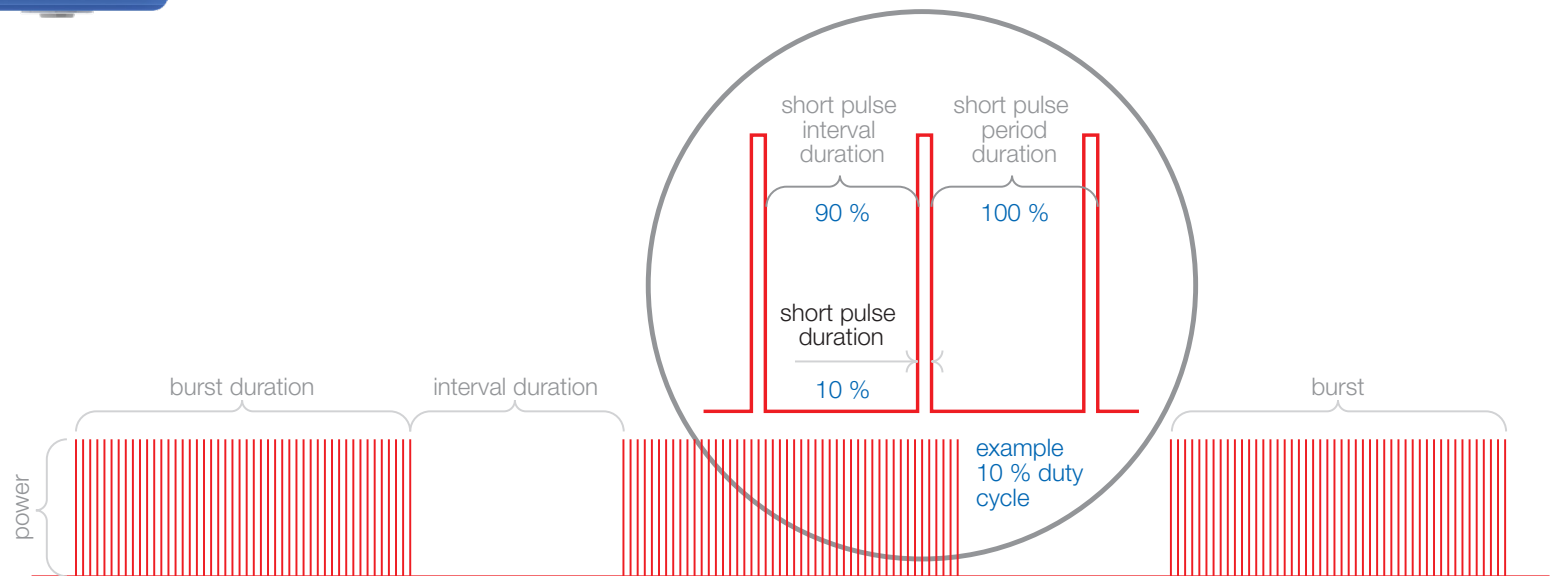


merilas
shortpulse **810**


STANDARD ACCESSORIES

- Foot switch
- Transport case

OPTIONAL ACCESSORIES

- Transscleral Cyclophotocoagulation probes for Glaucoma
- Other endoprobes
- Safety goggles

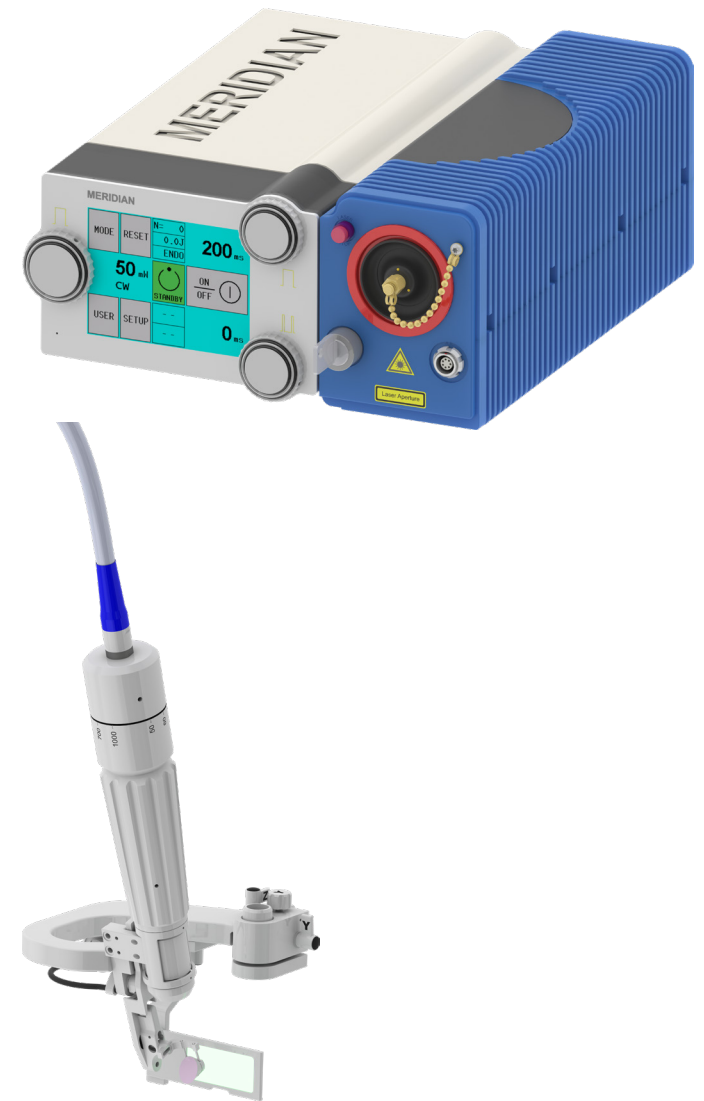

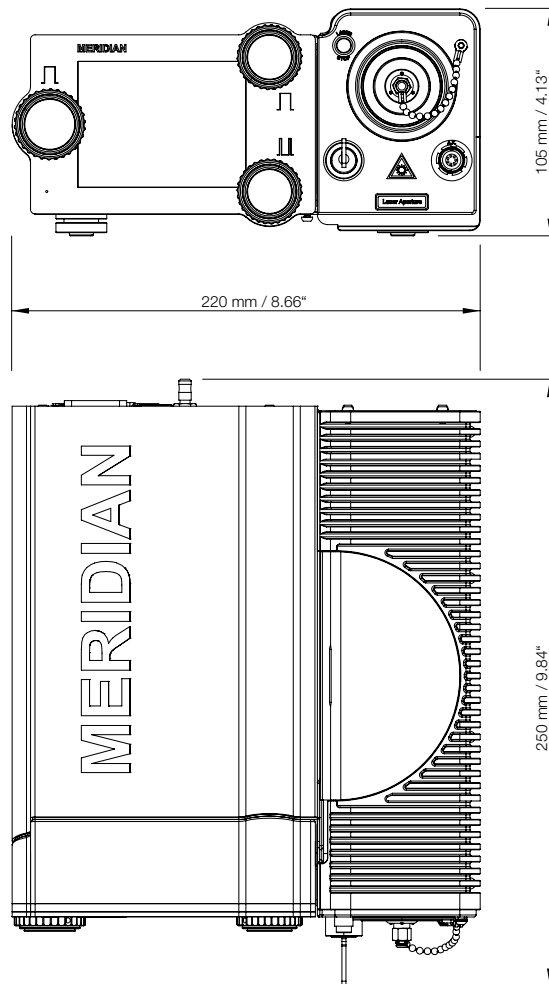


merilas

shortpulse **810**

TECHNICAL SPECIFICATIONS*

Device description	Merilas 810 shortpulse merilas shortpulse 810
Safety Classifications	Class 4
Wavelength	Infrared (810 nm)
Power Output	50 – 3000 mW
Pulse Duration	CW (continuous wave, chopped) 1 ms – 10 000 ms
Pulse Interval	1 ms – 10 000 ms
SP-Mode Settings	Shortpulse (continuous wave, chopped) Shortpulse duration: 0.01 – 9.5 ms Shortpulse interval: 0.1 – 9.5 ms
Cooling	TEC
Aiming Beam	Diode 635 nm, (0-1 mW in 9 Steps)
Dimensions	25.0 × 22.0 × 10.5 cm
Total Weight	7.0 kg
Power Requirements	100 – 240 V, 50/60 Hz, 2 A max.

* All technical specifications are subject to change without notice. In accordance with the international general safety standards: IEC 60601-1:2005/AMD1:2012, IEC 60601-1-2:014, MDD 93/42/EEC. The laser safety is in accordance with the international standards: IEC 60825-1:2014 and IEC 60601-2-22:2007/AMD1:2012.

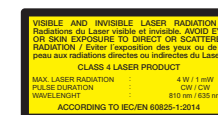

CE
0633

your laser specialist



Head office

Meridan Medical Group

Tel.: +41 33 334 11 11
Fax: +41 33 334 11 19
info@meridian.ch
www.meridian.ch

Switzerland

Meridian AG

Bierigutstrasse 7
CH-3608 Thun

Slovenia

Meridian Medical d.o.o

Plemljeva 8
1210 Ljubljana

Finland

Meridian Medical Oy

Elannotie 5
01510 Vantaa